QUICK START GUIDE

Read this manual carefully and save it for future reference.
Lisez attentivement ce manuel et conservez-le pour référence future.
Lea este manual con atención y guárdelo para futuras consultas.

IN THE BOX

Wyze Scale Ultra × 1
1.5V AAA Batteries × 4
Quick Start Guide × 1

DANS LA BOÎTE

Wyze Scale Ultra × 1
Piles AAA 1,5 V × 4
Guide de démarrage rapide × 1

EN LA CAJA

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Disclaimer: The products, accessories, user interfaces, etc. illustrated in this guide are intended for reference only. Due to product updates, the actual product and the illustrations in this guide may vary.

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
WELCOME TO WYZE BIENVENUE CHEZ WYZE BIENVENIDO A WYZE

Hi! Thanks for getting your hands on a Wyze Scale Ultra, and if this is your first Wyze product, welcome to the Wyze family! You just made a commitment to a better, more connected lifestyle. This guide will help you get things set up.

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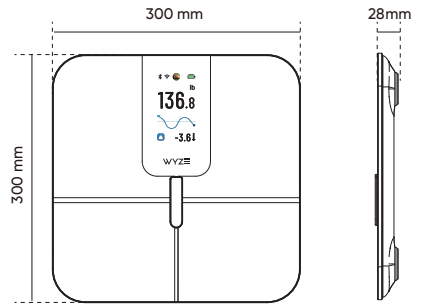SETUP CONFIGURATION CONFIGURACIÓN

➤ Before first use, open the battery compartment and insert the batteries.

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Note: When inserting the batteries, line up the positive and negative signs to make sure the batteries are installed correctly.

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➤ Download the Wyze app and sign in or create an account.

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
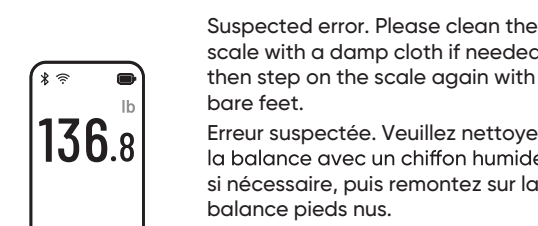
Suspected error. Please clean the scale with a damp cloth if needed, then step on the scale again with bare feet.

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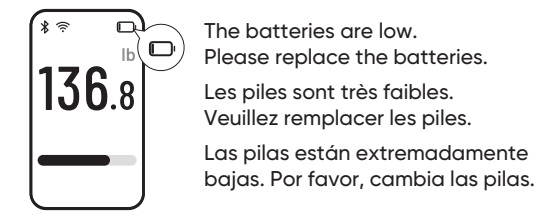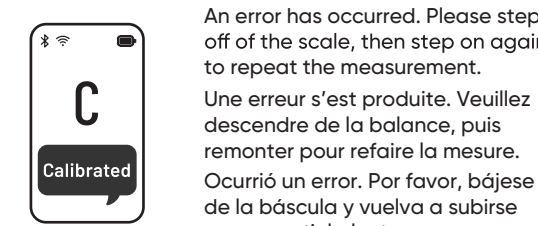
An error has occurred. Please step off of the scale, then step on again to repeat the measurement.

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Please recycle where facilities exist. Waste electrical products should not be disposed of with household waste. Check with your local authority or retailer for recycling advice.

Weight not in supported range. Please measure again. The poids n'est pas dans la fourchette prévue. Veuillez mesurer à nouveau. El peso está fuera del rango permitido. Por favor mide el peso otra vez.


EXPLANATION OF SYMBOLS ON UNIT EXPLICATION DES SYMBOLES SUR L'UNITÉ EXPLICACIÓN DE SÍMBOLOS EN LA UNIDAD

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➤ In the Wyze app, tap the + plus sign, then Device > Lifestyle > Wyze Scale Ultra. Follow the instructions in the Wyze app to complete setup. To connect, enable your phone's Bluetooth® and tap on the scale's surface to activate it.

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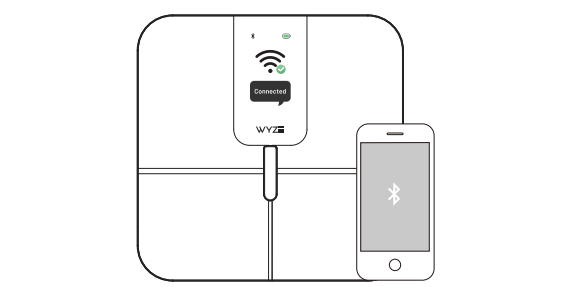➤ After the connection is established, place the scale on a hard, flat surface and step onto the scale with bare feet to begin the measurement. Make sure your feet are touching the 4 electrodes for the most accurate measurements.

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
CAUTION ATTENTION ATENCIÓN

People fitted with an implanted medical device such as a pacemaker, an electronic life support system such as artificial heart/lungs, or that use a portable electronic medical device such as an electrocardiograph, and women who are pregnant, must only use this product while wearing shoes, or in the "Only Measure Weight" setting. Please refer to the manufacturer's instructions for your medical device for more information. Do not use this product unless you are wearing shoes, or activate the "Only Measure Weight" setting, if you have an implanted medical device, have an electronic life support system, use a portable electrical medical device, or are pregnant. Pregnant women should also consult their healthcare provider before use.

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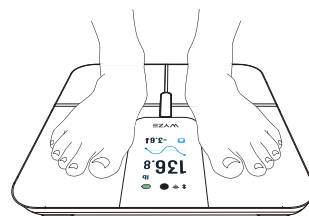On the back of your scale, find the UNIT button under the battery cover. Press it to change the units from lb to kg. To factory reset the scale, press and hold the same button for 5s until the TFT display shows "Clear".

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- Do not stand on the edge of scale, as you may fall over or influence the accuracy of the measurement.
- Do not use the scale on a tiled or wet floor, as this may result in a fall.
- When measuring your heart rate with your phone's camera lens, do not touch your phone's camera flash. Your phone's camera flash could overheat and burn your fingers.
- Make sure the surface of the scale is clean before you step onto it, to prevent any slippage.
- Treat your scale with care. Do not drop it or jump on it. The scale is designed to be stood on; misuse or abuse may break the electronic sensors, cause injury, or result in inaccurate measurements.

- Avoid using this scale near strong magnetic fields, such as microwave ovens, etc.
- This device complies with part 15 of the FCC Rules. Its operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.
- Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
- This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(S). Operation is subject to the following two conditions: (1) This device may not cause interference. (2) This device must accept any interference, including interference that may cause undesired operation of the device.
- This equipment complies with IC RSS-102 radiation exposure limits for residential uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.
- Bluetooth® associated logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Wyze Labs, Inc. is permitted under license. Other trademarks and trade names are those of their respective owners.

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NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against

harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: - Reorient or relocate the receiving antenna. - Increase the separation between the equipment and receiver. - Connect the equipment to an outlet on a circuit different from that to which the receiver is connected. - Consult the dealer or an experienced radio/TV technician for help.

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CARE & MAINTENANCE ENTretien Et Maintenance CUIDADO Y MANTENIMIENTO

- Avoid high temperatures and irradiation.
- If the scale is stored close to heating point, allow it to acclimate to room temperature before use.
- Remove the batteries if the scale is going to be used for 3 months or more.
- Clean the scale with a soft damp cloth. Do not use abrasive or solvent-based cleaners, and never immerse it in water.
- The scale will maintain its safety and performance features for at least 10,000 measurements.
- Évitez les températures élevées et l'irradiation.
- Si la balance est stockée près du point de congélation, laissez-la s'acclimater à la température ambiante

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FCC WARNING

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment. This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: - Reorient or relocate the receiving antenna. - Increase the separation between the equipment and receiver. - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected. - Consult the dealer or an experienced radio/TV technician for help.

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ISED / FCC RADIATION

EXPOSURE STATEMENT

This equipment complies with FCC/IC radiation exposure limits set forth for an uncontrolled environment. This equipment is installed and operated with minimum distance 14.23 mm between the radiator and human body.

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SPECIFICATIONS

Product Name: Wyze Scale Ultra Black
Model: WL_SCLUB
Product Size: 300 mm x 300 mm x 28 mm
Product Weight: 3.96 lbs / 1.80 kg
Power: 4 x 1.5V AAA batteries (included)
Display Units: lb/kg
Measurement Range: Weight: 11 - 400 lbs / 5 - 180 kg
Operating Temperature: 41°F - 104°F (5°C - 40°C)
Operating Humidity: 20 - 85% RH
Storage and Transport Temperature: -4°F - 131°F (-20°C - 55°C)

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