QUICK START GUIDE

Read this manual carefully and save it for future reference.
Lisez attentivement ce manuel et conservez-le pour référence future.
Lea este manual con atención y guárdalo para futuras consultas.

WELCOME TO WYZE BIENVENUE CHEZ WYZE BIENVENIDO A WYZE

Hi! Thanks for getting your hands on a Wyze Scale Ultra BodyScan, and if this is your first Wyze product, welcome to the Wyze family! You just made a commitment to a better, more connected lifestyle. This guide will help you get things set up.

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Disclaimer: The products, accessories, user interfaces, etc. illustrated in this guide are intended for reference only. Due to product updates, the actual product and the illustrations in this guide may vary.

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
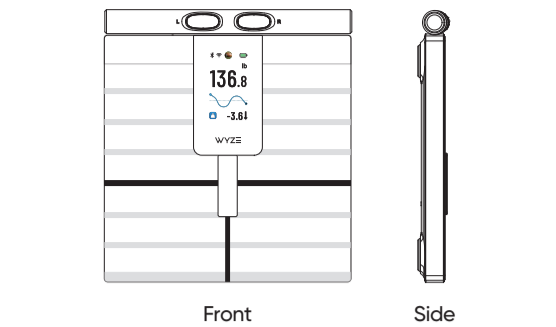
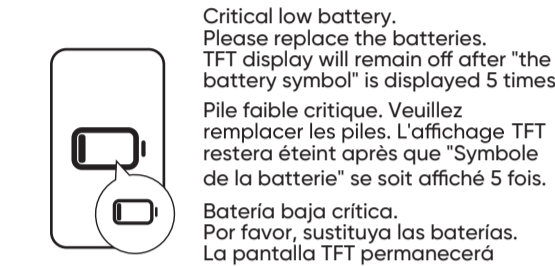
➔ In the Wyze app, tap the + plus sign. On the Add menu, tap **Device > Lifestyle > Wyze Scale Ultra BodyScan**. Follow the instructions in the Wyze app to complete setup. To connect, enable your phone's Bluetooth® and tap on the scale's surface to activate it.

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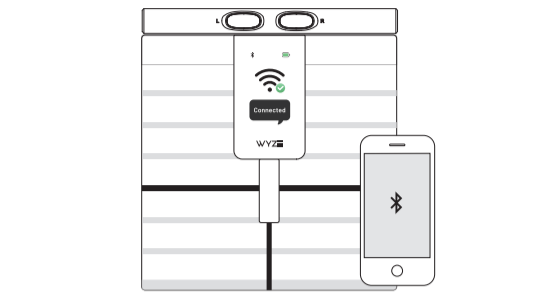
➔ After connecting to Bluetooth, step on the scale with bare feet and pull up the handle. Position your fingers on the electrodes and keep your arms straight out at a 30°-45° angle from your body.

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CAUTION ATTENTION ATENCIÓN

- People fitted with an implanted medical device such as a pacemaker, an electronic life support system such as artificial heart/lungs, or that use a portable electronic medical device such as an electrocardiograph, and women who are pregnant, must only use this product while wearing shoes, or in the "Only Measure Weight" setting. Please refer to the manufacturer's instructions for your medical device for more information. Do not use this product unless you are wearing shoes, or activate the "Only Measure Weight" setting, if you have an implanted medical device, have an electronic life support system, use a portable electrical medical device, or are pregnant. Pregnant women should also consult their healthcare provider before use.
- To activate the "Only Measure Weight" setting: 1. Open the Wyze App and tap Wyze Scale Ultra BodyScan. 2. Tap on the gear to open Scale Settings. 3. Tap "Users", then select your email address to open the "AccountInfo" screen. 4. Tap the slider next to "Only Measure Weight," then "Yes" to confirm.
- Some physical conditions could affect hydration levels that may lead to inaccurate results. Please consult your healthcare provider for more information.
- Only use 1.5V AAA batteries in the scale. Using the incorrect type of batteries creates a risk of explosion.
- Do not expose batteries to extremely high temperatures or fire. Doing so creates a risk of explosion or leakage of flammable liquid or gas.
- Do not manually crush, cut or deform batteries. Doing so creates a risk of explosion or leakage of flammable liquid or gas.
- Do not subject batteries to extremely low air pressure. Doing so creates a risk of explosion or leakage of flammable liquid or gas.
- Les personnes portant un dispositif médical implanté tel qu'un stimulateur cardiaque, un système de survie électronique tel qu'un cœur/poumons artificiels, ou qui utilisent un dispositif médical électronique portable tel qu'un électrocardiogramme, et les femmes enceintes, ne doivent utiliser ce produit qu'en portant chaussures ou dans le paramètre "Mesurer uniquement le poids".

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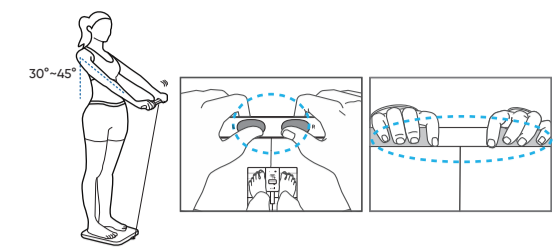On the back of your scale, find the UNIT button under the battery cover. Press it to change the units from lb to kg. To factory reset the scale, press and hold the same button for 5s until the TFT display shows "Clear".

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- Do not stand on the edge of scale, as you may fall over or influence the accuracy of the measurement.
- Do not use the scale on a tiled or wet floor, as this may result in a fall.
- When measuring your heart rate with your phone's camera lens, do not touch your phone's camera flash. Your phone's camera flash could overheat and burn your fingers.
- Make sure the surface of the scale is clean before you step onto it, to prevent any slippage.
- Treat your scale with care. Do not drop it or jump on it. The scale is designed to be stood on; misuse or abuse may break the electronic sensors, cause injury, or result in inaccurate measurements.

- Avoid using this scale near strong magnetic fields, such as microwave ovens, etc.
- This device complies with part 15 of the FCC Rules. Its operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.
- Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
- This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(S). Operation is subject to the following two conditions: (1)This device may not cause interference. (2)This device must accept any interference, including interference that may cause undesired operation of the device.
- This equipment complies with IC RSS-102 radiation exposure limits set forth in an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.
- Bluetooth® associated logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Wyze Labs, Inc. is permitted under certain other trademarks and trade names are those of their respective owners.

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harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: - Reorient or relocate the receiving antenna. - Increase the separation between the equipment and receiver. - Connect the equipment to an outlet on a circuit different from that to which the receiver is connected. - Consult the dealer or an experienced radio/TV technician for help.

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CARE & MAINTENANCE ENTRETIEN ET MAINTIENANCE CUIDADO Y MANTENIMIENTO

- Avoid high temperatures and irradiation.
- If the scale is stored close to freezing point, allow it to acclimate to room temperature before use.
- Remove the batteries if the scale is not going to be used for 3 months or more.
- Clean the scale with a soft damp cloth. Do not use abrasive or solvent-based cleaners, and never immerse it in water.
- The scale will maintain its safety and performance features for at least 10,000 measurements.
- Évitez les températures élevées et l'irradiation.
- Si la balance est stockée près du point de congélation, laissez-la s'acclimater à la température ambiante

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FCC Compliance Statements

Caution: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: -- Reorient or relocate the receiving antenna. -- Increase the separation between the equipment and receiver. -- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected. -- Consult the dealer or an experienced radio/TV technician for help.

ISED Compliance Statements

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(S). Operation is subject to the following two conditions: (1) This device may not cause interference. (2) This device must accept any interference, including interference that may cause undesired operation of the device.

RF Exposure Information (SAR)

Specific Absorption Rate (SAR) refers to the rate at which the body absorbs RF energy. This equipment has been tested and meets applicable limits for radio frequency (RF) exposure. During testing, device racks are set to their highest transmission levels and placed in positions that simulate use near the body or limb, with 0 mm separation. Accessories with metal parts may change the RF performance of the device and use such accessories should be avoided. The SAR limit is 4.0 W/kg over 10 gram of tissue in US/Canada.

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SPECIFICATIONS

Product Name: Wyze Scale Ultra BodyScan
Model: WL_SCLETB/WL_SCLETW
Product Size: 13.39 x 15.43 x 2.13 in / 340 x 392 x 54 mm
Product Weight: 5.95 lbs / 2.7 kg
Power: 4 x 1.5V AAA batteries (included)
Display: 1.5 in / 3.8 cm
Measurement Range: Weight: 11 - 400 lbs / 5 - 180 kg
Operating Temperature: 41°F - 104°F (5°C - 40°C)
Operating Humidity: 20 - 85% RH
Storage and Transport Humidity: ≤90% RH
Storage and Transport Temperature: -4°F - 131°F (-20°C - 55°C)

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Apple logo and App Store are trademarks of Apple Inc., registered in the U.S. and other countries and regions. Google Assistant, Google Home, Google Play, and the Google Play logo are trademarks of Google LLC.

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Wyze and Wyze Scale Ultra BodyScan are trademarks of Wyze Labs, Inc.
Wyze et Wyze Scale Ultra BodyScan sont des marques déposées de Wyze Labs, Inc.
Wyze y Wyze Scale Ultra BodyScan son marcas comerciales de Wyze Labs, Inc.